

Procenta[®]

Placental-Derived Allograft



Procenta[®] is an acellular placental-derived graft made from soft connective tissue that is hydrophilic, sterile and requires no re-hydration or orientation.



PROCENTA'S CAPABILITIES

- Provides a conformable, adherent scaffold to cover a site of injury
- Protects the wound site from surrounding environment
- Retains fluid at the wound site



PRODUCT CHARACTERISTICS

- Procenta[®] is an acellular, highly hydrophilic, structural extracellular matrix derived from placental connective tissue.
- Structural, human tissue allograft (FDA/TRG recommendation letter received).
- Comprised of extracellular matrix scaffold that serves as a cover, to offer protection from the surrounding environment, and to retain fluid.
- Allograft produced from donated human placental tissue. Intended for homologous use to supplement or replace damaged or inadequate connective tissue.
- Stored at ambient temperature, ready for application in a sterile vial.
- 4 year shelf life.


Disclaimer:

Procenta[®] is not intended to diagnose, treat, cure or prevent any disease, and should not be used for treatment of COVID-19.

CONTACT [Eric Karukin call or text 727-645-7774](tel:727-645-7774)
eric@products4doctors.com



PRODUCT PACKAGING & CONTENT

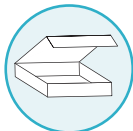

Outer Box that contains sterile packaged Procenta® vial.



Vial is packaged in a single pouch.



Allograft Tracking Card that must be completed and returned to Russell Health.



Allograft Return Card must be completed and sent back to Russell Health in case product is returned.

PRODUCT APPLICATION



Open product package & remove the Tyvek peel pouch containing the vial.



Using aseptic technique, peel open the Tyvek peel pouch and present sterile vial onto sterile field.



Open the vial and apply product as needed.

PRODUCT PORTFOLIO

Catalog #	Product Description
PRO100	Procenta® Placental-Derived Allograft 100mg
PRO200	Procenta® Placental-Derived Allograft 200mg
PRO300	Procenta® Placental-Derived Allograft 300mg



AMBIENT VERSION



DIRECTLY APPLIED



4 YEAR SHELF LIFE

PROCESSING

The manufacturing of Procenta® Tissue-based Allograft does not involve the combination of the placental tissue with another article, except for a placental storage agent.

After minimally manipulated processing, Procenta® is terminally sterilized using E-beam irradiation to ensure recipient safety.

The Procenta® Allograft processor is Lucina BioSciences, LLC. A proprietary, minimally manipulative manufacturing process (Nativus) ensures the original relevant characteristics of the tissue are retained.

The tissue bank is registered with the Food and Drug Administration (FDA).



Disclaimer:

Procenta® is not intended to diagnose, treat, cure or prevent any disease, and should not be used for treatment of COVID-19.